


ECHELON™ 3000 Stapler

Optimized device performance guide


Loading the device

Proper reload selection: Tissue thickness should be carefully evaluated prior to using the instrument. Refer to the **Reload Product Codes Chart**, located in the Instructions for Use, for proper reload selection. Selection of the appropriate staple reload should be based upon the combined thickness of both the tissue and the staple line reinforcement material if used.

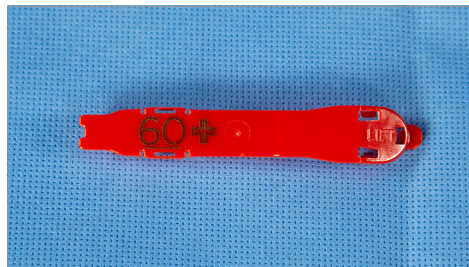



Step 1

Verify that the reload size matches the instrument size being used (e.g., use 60mm reload with 60mm stapler).

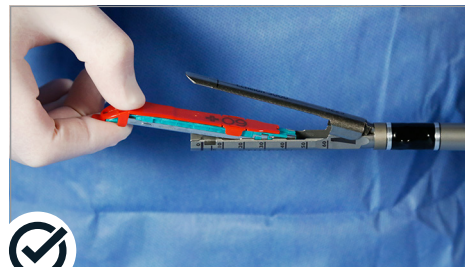


Step 2

Ensure the staple retaining cap is present during loading to minimize the potential for premature sled movement, which could lead to inadvertent lockout.

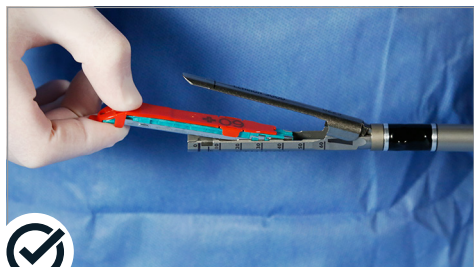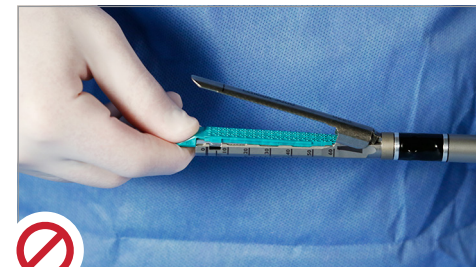


Step 3

Insert the reload by sliding it against the bottom of the reload jaw until the stapler alignment tab stops in the reload alignment slot of the jaw. Positioning the reload at the bottom of the reload jaw could prevent inadvertent sled movement and lockout.

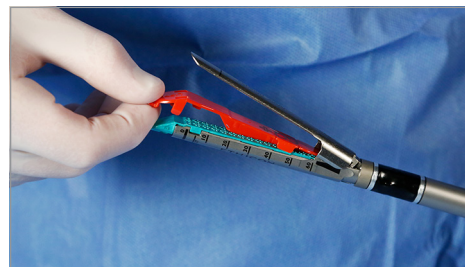

Step 4

Remove the staple retaining cap by lifting it straight up and discard.

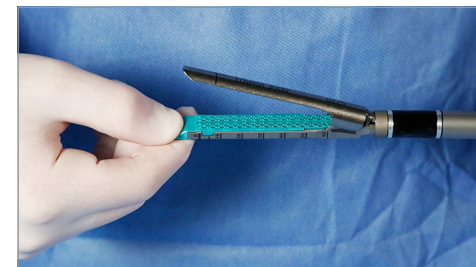


Step 5

Observe the surface of the new reload to ensure no drivers are visible and that the reload lies flat in the jaw channel.

ECHELON™ 3000 Stapler


Optimized device performance guide

Device articulation and centering

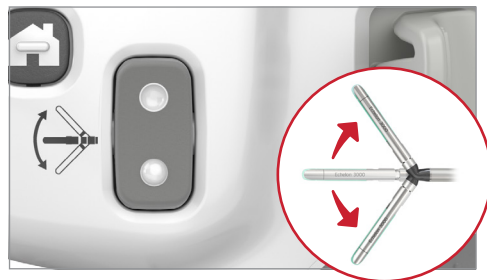

Step 1

After inserting the battery, press/test the left and right articulation buttons (buttons w/ white circles) to confirm powered articulation is ready for use.

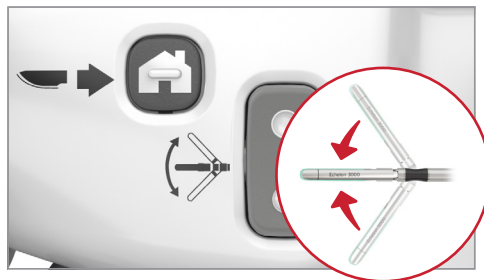

Step 2

Press the home button to return the jaws to the home (straight) position.




Step 3

Practice with the Anvil Jaw up and the Anvil Jaw down. The device senses the orientation of the jaw relative to the handle to keep the articulation direction controls consistent. Once the jaws have reached the maximum angle in either direction, if the Articulation Control Button is pressed again, the instrument will provide audible/haptic feedback to indicate that the jaws cannot be articulated further in that direction.




Partial close articulation

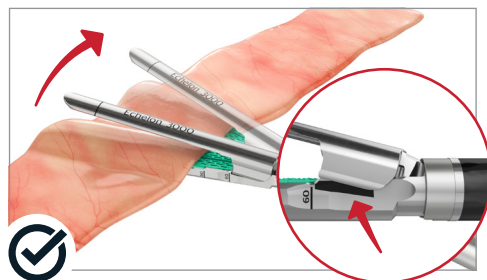

Step 1

The device allows for articulation to occur when the jaws are open and the black partial close indicator (denoted by the red arrow) is visible.

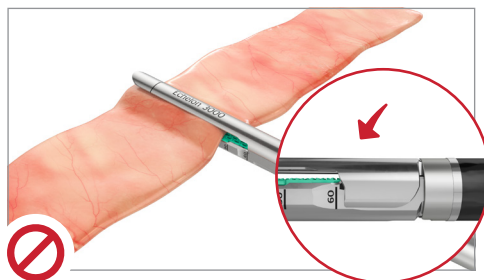

Step 2

When closing the jaws, once the base of the Anvil Jaw covers up the black partial close indicator on the side of the reload jaw, the device is no longer able to be articulated.




If the articulation buttons are pressed when the jaws are closed, the device will provide haptic feedback to the user to indicate that it cannot be articulated while the jaws are closed.

ECHELON™ 3000 Stapler

Optimized device performance guide

Device closure

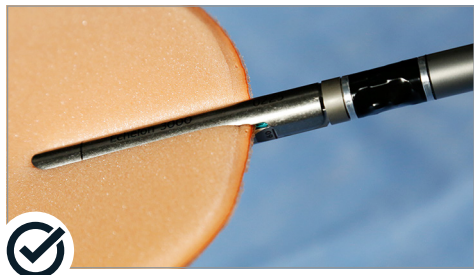


Note: The user may experience audible / haptic feedback during clamping and unclamping of the device. This is typical of device operation.

Step 1


Position the instrument around the tissue to be stapled. When closing the device, ensure that the tissue lies flat and is positioned properly between the jaws. Any “bunching” of tissue along the reload, particularly in the crotch of the jaws, may result in an incomplete staple line or inadvertent lockout. Holding the jaws in place for 15 seconds after closing and prior to firing may result in better compression and staple formation.

Firing the device

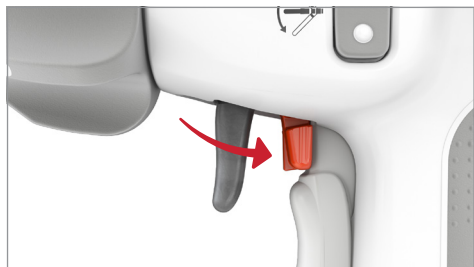


Step 1

Pull back the Red Firing Trigger Lock to enable the Dark Gray Firing Trigger to be pulled.

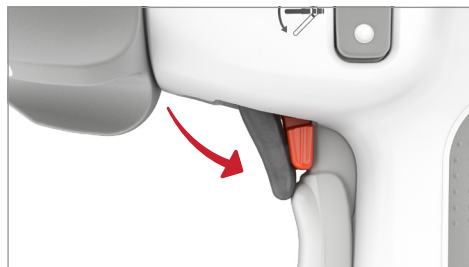

Step 2

Fire the instrument by pulling the firing trigger. Continue to depress the trigger until the motor stops. To complete the firing sequence, release the firing trigger to automatically return the knife to the home position.

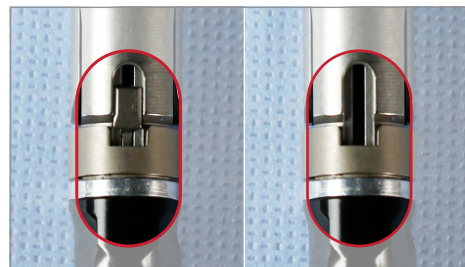

Knife Blade
Visible

Knife Blade
Not Visible


Note: Alternatively, for enhanced tissue compression when using the ECHELON™ 3000 Stapler, the user may pulse the firing trigger throughout the firing cycle until the transection is completed.

Similarly, if the firing sequence is interrupted, the user may pulse the firing trigger to continue the transection. It is recommended that the pulsing technique not be initiated until after the knife blade indicator has moved into the closed bottom channel and is no longer visible to minimize potential of inadvertent lockout.

ECHELON™ 3000 Stapler


Optimized device performance guide

Opening the jaws

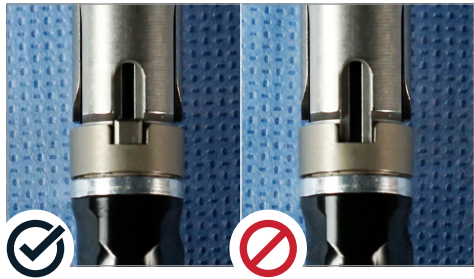


Step 1

Ensure the knife is in the home position by observing the knife blade indicator under the cartridge jaw.

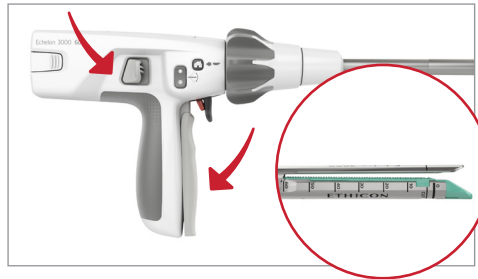

Step 2

To open the jaws, first squeeze the closing trigger, and then, without releasing the closing trigger, push the clamp release button on the side of the device. Maintaining pressure on the clamp release button, slowly release the closing trigger. If the jaws do not automatically open after the clamp release button is pressed and the knife is in the home position, then gently pull the closing trigger upward (away from the handle) until both firing and closing triggers return to their original positions.

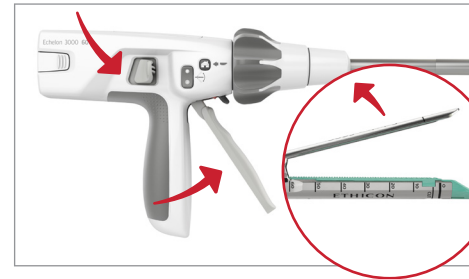

Caution: For insertion and removal of articulating instruments, the jaws of the instrument must be straight, parallel to the shaft of the instrument. Failure to have the instrument jaws in the straight position will result in difficult insertion or withdrawal of the instrument and may result in damage to the instrument.

Preparation for additional firings



Step 1

Prior to reloading the stapler, hold the instrument in a vertical position with anvil and cartridge jaw completely submerged in sterile solution. Swish vigorously to remove any loose formed staples.



Step 2

Wipe the inside and outside surfaces of the anvil and cartridge jaws to clean any unused staples from the instrument.

ECHELON™ 3000 Stapler

Optimized device performance guide


Jaws will not open

If the instrument is fired completely, the knife will return to the home position automatically. At any time, if the knife does not return to the home position and the jaws will not open:



Step 1

Home button: Ensure the Battery Pack is securely installed and the instrument has power; then, press the home button on either side of the device to return the knife to the home position.

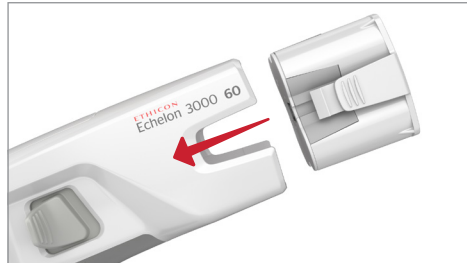

Step 2

Battery: If pressing the home button is not successful, remove and re-install the battery to check that it is securely installed and that the device has power.
Repeat step 1: Press the home button to return the knife to the home position.



Step 3a

Manual override: If the knife does not return, use the manual override. First, remove the access panel labeled "Manual Override" on the top of the handle.



Step 3b

Move the lever FORWARD AND BACKWARD until it can no longer be moved. The forward and backward movement of the lever may require multiple cycles and should be completed until the knife blade is confirmed to be in the home position.

Note: After the manual override is used, the device becomes disabled for any subsequent firings.